


KEY CONCEPT: FLEXIBLE LATERAL/MEDIAL FOREFOOT CONTAINMENT,  
CONTROLS MILD TO MODERATE  
PRONATION/SUPINATION

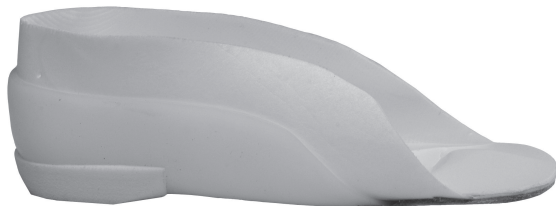

The DAFO 5 Softy uses a soft foam liner to provide comfort and flexible containment of the lateral/medial forefoot. The lower plastic trimline provides good hindfoot stability without compromising comfort. The heel remains stable in the corrected position while the soft foam liner assists in containing a mild to moderately pronated or supinated foot. The plantar surface flexible plastic shell extends distal to the Met Heads, transitioning into a flexible leather sole. The soft foam liner is ideal for patients with challenging foot anatomy, difficult to correct feet, or those prone to pressure points and irritation.

**RECOMMENDED FOR PATIENTS WHO:**

1. Exhibit mild to moderate pronation or supination with adequate ankle control
2. Have mild heel eversion with forefoot abduction
3. Are larger/more mature
4. Have difficult to correct feet or feet with sensitive bony prominences.

This DAFO 5 Softy is shown with  
a soft foam inner liner.

**Hindfoot:** Stable

**Forefoot:** Medial/Lateral Containment

**Ankle:** Plantarflexion Free/Dorsiflexion Free